## **Genome Browser**

## Table of Contents:

- Genome browser is a genomic viewer that allows to visuazile in different zones of visualization(trakcs) different features (genes, snp, variations...).
- Architecture
- Tracks and Panels

## Architecture

The genome browser is composed of 2 panels showing the careotype and chromosomes of the selected species (blue) and in the lower part different tracks or zones of visualization of genomic features (green).

0 9-	Homos	agiens •	13+	h	×≈	M	-	_		•	+ M	6 1	17		33.32	994308	32999	64	Get	×		• 2	r (0			۹									
1	1				5								53				1				17			ļ			-		22	III.	ļ		17	ņ	
1	-		ą		100	4111	100	1 8	401	432	(1))	1110	040	1	042	2			1	61136	0110	6110	100		ą	- i		612	- 10	1	8			600	ŝ.
egion over	view															Wr		te: 117	5Hs																
OPERATION Concernent																	32	963MI																	329940
Sectioned	sew -	**													-		_			-															
etailed inf	omatio															W		be: 140 26381	Tats.																32776
	AACA		ATTGA	CAAC		ААТС	TAG	стте	стта		IGA <b>C</b> A	ATT	SAA7	тс	AGG	TTAA	TTA		стт	ATA	100 <b>0</b>	TCTI	677	****	CAT	ГАТА	сало	GAT	слаа	ACCAN	астс	ITA		TTCAC	CACATO
iene =	~~																_																		
Variation	- • •				- 11415 1555555	200 r A		n7 - A				18365	1726.A	1	614021	122 C > A		10000	274067 02956 ( 129633	1 × 1					902734	077×0	nan g Pi	472201	1>T 84.C+T	-	мыс-				11112ACEA 1844711735- 175634255

The panels are defined when you instantiate the Genome Browser, the tracks are added later.

One of the main characteristics is their joint movement. When you move to a zone of the genome all the tracks move together to visualize the same zone of genome.

The navigation bar is shown at the top. Ih this bar you can select settings that will affect the entire genome browser, such as the species selector. From left to right it shows:

- Recharge button: This button reloads the genome browser locating it in the initial position.
- Region history button: In this button you can reposition your genome browser in a position previosly searched.
- Species selector: In this selector the user can change species. The available species are those available in the cellbase database. (Cellbase species)
- Chromosome selector: This selector is a quick way to change the chromosome.
- Toggle Karyotype panel: Show or hide the karyotype panel.
- Toggle chromosome panel: Show or hide the chromosome panel.
- Toggle overview panel: Show or hide the overview panel.
- Min button: Minimun windows size.
- "-" button: In this button the user can decrease the zoom by making the area to be displayed larger.
- Zoom scrollbar: Scrollbar to adjusted to the desired zoom by the user.
- "+" button: In this button the user can increase the zoom by making the area appear smaller.
- Max button: Maximun windows size.
- · Windows size button: Show the window size in nucleotides.
- Position: This display shows the region where the tracks are located. The structure of region is "chromosome:start-end". This input is editable, the user can write the desired position and then click on the go button to change the region.
- Gol: This button positions the genome browser in the position shown in the display position.
- <<, <, >, >>: These buttons are to scroll through the genome browser. If the user wishes instead
  of clicking and dragging on the track to move, you can press these buttons to move laaterally
  within the genome. The "<<" and ">>" buttons will make a big jump in their directions while the
  movement when pressing "<" or ">" will be smaller.
- Compact button: if this button is pressed the tracks will be shown in extended mode occupying all the space necessary to show the information contined in its entirety. If the botton is not pressed, the tracks will have the predefined height in their creation and a lateral scroll bar will be shown when the information they contain is greater than the size they have.
- Gene display and search: For more comfort the user genome browser can be located in the region that occupies a certain gene. For this, the user can enter a gene in the display and give the magnifying glass to search.

## **Tracks and Panels**


• Karyotype: All human chromosomes are shown. The current display area is indicated by an orange line.

	Homo	uapierus *	11× 1	t X To	Min =	_	-	+ Max	147	12	22996302-	2996448	Gal	« < >	Quine .	Q							
Karyotype																							
Ţ				1										13	ļ		ļ	1	U.	ų	se.	ţ.	

 Chromosome: The selected chromosome is displayed. The current display area is indicated by an orange line.



• Gene overview: Region genes shown. This area is less than the area zoom genes. It is to have a global vision to the user that there are genes around.

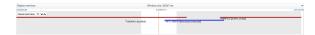

• Sequence: Bases of chromosomal sequence. If the viewing area is very large sequence it will not be shown.

Detailed information	Window size: 137 nts	
Sequence - VA	1004608	166457
GGCACCTTCCCCCATCCGTAGCCCCTCCCCAAGCGCGCTT	GTEEGEAGGGETGGAETGEGAEAEEGAEGTEEEGAAGAEGATEEAGETGGTTTGGEGGEAG	CACTGGGGAATGGTGCAGACGGAGGCTCCGTACA/

• Gene: Gene and transcripts found in that region.

Gene – VA	
	00 00 00 jamanajandikrifiziki

In a closer view we can see exons painted red when they are protein coding. When you place the mouse on top of the exon, the track show the amino acid and a tooltip which exon description.

o Pro Pro Arg Pro Pro Arg Cin Arg Cin Gin Gin Gin Gin Pro Cin Leu Ser Leu Arg Giy P		n - ENSERVALISATING	
	Trans Hit dei das her	onto: EnvST00000800200 Aurpool: unceffered exception: unceffered envol: 320744864-30002315 {} rgth: 87,465	
	inn dry pro- edit coli coli coli coli coli coli coli coli	yeb, 276 metric fill and the second	

• SNP: SNP found in that region.

SNP - VA	
7109753402 70271314394	
10513403	

If the size of the region is no longer, the snp be displayed , if not a histogram of these .

SNP 🖬 = 🗸 🗛	

• Variant browser - compact mode: In this mode the variants are displayed in mode browser.

Veriet - v A						
10-A	12.229	1659270	ra	2/2/10/06110/2	12.22000581 5-0	i.
	13329983624.0		13.32998891.01			i.

• Variant genotype - extended mode: In this mode there are samples. The genotypes of the samples are drawn, showing in blue 0/0, in yellow 0/1 and red 1/1.

• Bam: The coverage is shown at the top. In the lower part the readings are shown.

Algenet - VA X